

V3 SECTIONAL MATRIX SYSTEM

KEY BENEFITS & STARTER PACK CONFIGURATION



PERFECT CONTACTS - IDEAL ANATOMY

- + Super-elastic nickel-titanium spring for optimal separation force and tight contacts
- + Wide indications for use - glass-fiber tines will not collapse into wide cavities, even when a cusp is missing
- + Tines anatomically shaped for superior matrix retention and adaptation - and minimal finishing
- + Tines' V-shape accommodates the wedge
- + V3 Tab-Matrix with pronounced marginal ridge and vertical profile for correct anatomical form
- + Wave-Wedge stackable, self-locating and self-seating for optimal seal on matrix
- + Spring tilt allows easy stacking of multiple rings
- + Ring very stable on undercuts, and in forceps during placement

V3 STARTER PACK INCLUDES:

- + 2 x V3 Ring Universal (green)
- + 1 x V3 Ring Narrow (yellow)
- + 25 x Wave-Wedges Small (white)
- + 25 x Wave-Wedges Medium (pink)
- + 25 x Wave-Wedges Large (purple)
- + 25 x V3 Tab-Matrices 4.5mm
- + 25 x V3 Tab-Matrices 5.5mm
- + 25 x V3 Tab-Matrices 6.5mm
- + 5 x WedgeGuard Small (white)
- + 5 x WedgeGuard Medium (pink)
- + 5 x WedgeGuard Large (purple)
- + 1 x Pin-Tweezers
- + 1 x Forceps

PLUS INCLUDES SUPERCURVE NON-STICK MATRIX SAMPLES:

- + 10 x SuperCurve Non-Stick Matrices 4.5mm
- + 25 x SuperCurve Non-Stick Matrices 5.5mm
- + 25 x SuperCurve Non-Stick Matrices 6.5mm

PRODUCT CODE: 403318



www.trident.com
supercurve
 NON-STICK matrix

trident
 Innovative. Simple. Smart.

Visit www.suredental.com for technique video



Call for pricing


- A) New shape with pronounced marginal ridge.
- B) An S- shaped vertical contour provides close replication of the true anatomy of the proximal surface and allows the matrix to slide past the gingival margin.
- C) V3 matrices are 30° thin
- D) Greater curvature on the horizontal plane means more of the curvature remains as the matrix wraps around the tooth.
- E) The Tab acts as a handle for easy placement.
- F) Holes in the tabs of the V3 matrix allow easy removal with pin tweezers.

The V3 Tab-Matrix takes the Tab-Matrix a step further, introducing a pronounced ridge on the occlusal margin. This curved edge produces completely natural-looking Class II composite restorations that require minimal finishing, saving time, money and effort.

Available in pkgs 50 or 100 4.5mm,5.5mm,6.5mm

SuperCurve

"The new SuperCurve matrix from Triodent is the best matrix I have ever used. Its shape makes it incredibly easy to place, and creates perfect contours with minimal finishing. The non-stick surface makes removal a breeze. Finally, a perfect matrix!"

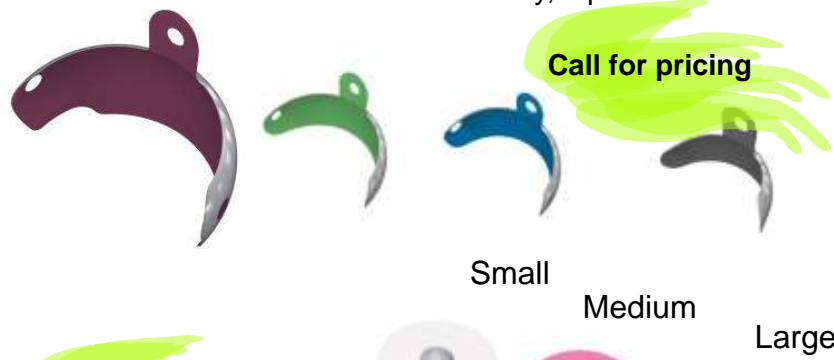

Call for pricing

Small Medium Large

Original V-Matrix

The Tab-Matrix is anatomically accurate and easy to handle, thanks to the tab and holes which fit Triodent's Pin-Tweezers, making placement and removal simple



\$75.00
100 Pk

Wave-Wedge

Used with the V3 Ring, the Wave-Wedge has unsurpassed sealing capability at the gingival margin. Its smart features, including a central concavity - the 'wave' - ensure it stands head and shoulders above any other wedge available.

Available in pkgs 100



Call for pricing



SureDental™

1-888-277-2335